

Evaluation and Comparison of FTA Card and CTAB DNA Extraction Methods for Non-Agricultural Taxa

Authors: Siegel, Chloe S., Stevenson, Florence O., and Zimmer,

Elizabeth A.

Source: Applications in Plant Sciences, 5(2)

Published By: Botanical Society of America

URL: https://doi.org/10.3732/apps.1600109

BioOne Complete (complete.BioOne.org) is a full-text database of 200 subscribed and open-access titles in the biological, ecological, and environmental sciences published by nonprofit societies, associations, museums, institutions, and presses.

Your use of this PDF, the BioOne Complete website, and all posted and associated content indicates your acceptance of BioOne's Terms of Use, available at www.bioone.org/terms-of-use.

Usage of BioOne Complete content is strictly limited to personal, educational, and non - commercial use. Commercial inquiries or rights and permissions requests should be directed to the individual publisher as copyright holder.

BioOne sees sustainable scholarly publishing as an inherently collaborative enterprise connecting authors, nonprofit publishers, academic institutions, research libraries, and research funders in the common goal of maximizing access to critical research.

APPLICATION ARTICLE

EVALUATION AND COMPARISON OF FTA CARD AND CTAB DNA EXTRACTION METHODS FOR NON-AGRICULTURAL TAXA¹

CHLOE S. SIEGEL^{2,5}, FLORENCE O. STEVENSON^{3,5}, AND ELIZABETH A. ZIMMER^{4,6}

²Department of Crop Sciences: Plant Biotechnology and Molecular Biology, University of Illinois at Urbana-Champaign, 1102 S. Goodwin Avenue, Urbana, Illinois 61801 USA; ³Department of Mechanical Engineering, Department of Mathematics, University of Maryland, 2181 Glenn L. Martin Hall, Building 088, College Park, Maryland 20742 USA; and ⁴Department of Botany, Smithsonian Institution National Museum of Natural History, NMNH-MRC, 166 Smithsonian Institution, P.O. Box 37012, Washington, D.C. 20013-7012 USA

- Premise of the study: An efficient, effective DNA extraction method is necessary for comprehensive analysis of plant genomes. This study analyzed the quality of DNA obtained using paper FTA cards prepared directly in the field when compared to the more traditional cetyltrimethylammonium bromide (CTAB)—based extraction methods from silica-dried samples.
- Methods: DNA was extracted using FTA cards according to the manufacturer's protocol. In parallel, CTAB-based extractions
 were done using the automated AutoGen DNA isolation system. DNA quality for both methods was determined for 15
 non-agricultural species collected in situ, by gel separation, spectrophotometry, fluorometry, and successful amplification and
 sequencing of nuclear and chloroplast gene markers.
- Results: The FTA card extraction method yielded less concentrated, but also less fragmented samples than the CTAB-based technique. The card-extracted samples provided DNA that could be successfully amplified and sequenced. The FTA cards are also useful because the collected samples do not require refrigeration, extensive laboratory expertise, or as many hazardous chemicals as extractions using the CTAB-based technique.
- *Discussion:* The relative success of the FTA card method in our study suggested that this method could be a valuable tool for studies in plant population genetics and conservation biology that may involve screening of hundreds of individual plants. The FTA cards, like the silica gel samples, do not contain plant material capable of propagation, and therefore do not require permits from the U.S. Department of Agriculture (USDA) Animal and Plant Health Inspection Service (APHIS) for transportation.

Key words: CTAB-based technique; DNA extraction, amplification, and sequencing; fluorometry; FTA cards; gel electrophoresis; spectrophotometry.

Successful extraction of DNA from plant tissues generally has been most successful from freshly field-collected tissue placed into either liquid nitrogen or, more commonly, into silica gel desiccant. When possible, the samples are then stored frozen in ultracold liquid nitrogen tanks or in boxes containing additional silica gel at -20°C to -80°C. Subsequent nucleic acid extractions most commonly have involved either some form of a cetyltrimethylammonium bromide (CTAB)-based technique (Doyle and Doyle, 1987) or use of a QIAGEN DNeasy kit (QIAGEN, Venlo, The Netherlands). However, previous research has also shown success with an alternative genomic DNA extraction method using Whatman FTA PlantSaver Cards

¹Manuscript received 16 September 2016; revision accepted 10 January 2017.

The study could not have been completed without the field expertise of Dr. Carol Kelloff. The authors would also like to thank Gabriel Johnson of the Smithsonian Institution and Juannan Zhou of the University of Maryland, who provided invaluable insights and mentorship throughout the research process. Laboratory and computer work were conducted in and with the support of the Laboratories of Analytic Biology facilities of the Smithsonian National Museum of Natural History.

⁵These authors contributed equally to this work.

⁶Author for correspondence: zimmerl@si.edu

doi:10.3732/apps.1600109

(Whatman, Maidstone, United Kingdom), with effectiveness for agricultural plant species, entomology, mycology, and other food sciences (Suzuki et al., 2006; Marques et al., 2010; Adugna et al., 2011; Bujang et al., 2011; Chandrashekara et al., 2012). Nevertheless, remarkably little research has been done to test the effectiveness of this extraction method on non-agricultural plant species. Additionally, it has been shown that methods employing DNeasy kits and/or CTAB-based methods may not be optimal for all plant genera and families because inhibitors can be coprecipitated with the DNA (Bustin and Nolan, 2004; Adugna et al., 2011). Therefore, it is proposed that alternative extraction methods, such as use of the FTA cards, be tested as to whether they are more effective for some plant genera.

To help gain a better understanding of the potential uses and limitations of FTA card extraction, this study assessed FTA card–extracted DNA quality in terms of concentration, spectral absorption, degree of fragmentation, and amplification and sequencing ability on a wide phylogenetic range of non-agricultural species. The data collected using these methods helped to indicate which plant species may be most compatible with the FTA card extraction method. Both successful and failed extractions provided valuable insights into the potential advantages and limitations of this alternative extraction method.

Applications in Plant Sciences 2017 5(2): 1600109; http://www.bioone.org/loi/apps © 2017 Siegel et al. Published by the Botanical Society of America. This open access article is a U.S. Government work distributed under the terms of the Creative Commons Public Domain Dedication (CC0 1.0).

MATERIALS AND METHODS

Collection—Samples from 15 phylogenetically diverse taxa possessing varying leaf characteristics were collected from the following vascular plant families: Apocynaceae, Aquifoliaceae, Asteraceae, Cactaceae, Cyperaceae, Fabaceae, Lamiaceae, Magnoliaceae, Oleaceae, Oxalidaceae, Poaceae, Typhaceae, Vitaceae, Pinaceae, and Aspleniaceae (Table 1). Collection of samples was completed on the morning of 23 June 2015, alongside a gravel utility road, sloping hillsides, and free-standing trees in Suitland, Maryland, USA. Voucher specimens for each species were deposited in the National Museum of Natural History's Herbarium (US), the University of Illinois at Urbana-Champaign Plant Biology Herbarium (ILL), and the Chicago Botanic Garden (CHIC).

Samples were preserved using Whatman FTA PlantSaver Cards and in silica gel (Flower Drying Type A silica with indicator; AGM Container Controls, Tucson, Arizona, USA). Samples were applied to the FTA cards directly in the field, then stored at room temperature. To make the plant print, a ceramic pestle was used like a hammer to smash the leaf tissue onto the card paper (Fig. 1). Enough plant prints were made for seven replicate extractions to be performed for each species.

Additionally, 2-3 in² of leaf tissue were collected into individual silica gelcontaining bags to be extracted later, using the CTAB-based technique on the AutoGen DNA isolation system (AutoGen, Holliston, Massachusetts, USA). Bags were stored in a plastic box containing additional silica gel at room temperature prior to extraction and then moved to a -20° C freezer for long-term storage.

Extraction

FTA card extraction—After collecting the samples on the FTA cards, small punched disks were removed from the sample cards followed by a series of washes on the disks. The protocol used was closely modeled after the method outlined in Adugna et al. (2011), with a modified disk size and centrifugation added after the incubation period to reduce bubbles. Eight disks, 2.0 mm in diameter, were used in this study so as to have a comparable total disk surface area to that found in Adugna et al. (2011). Additionally, the protocol of Adugna et al. called for centrifugation only after the addition of TE in the last step of the protocol. In our study, plates were centrifuged at 6200 rpm for 2 min at 4°C in between each purification wash.

Using the 2.0-mm-diameter Harris Micro Punch and Mat (Whatman), 56 disks were removed from each specimen's plant-pressed FTA card. A 96-well, Square V-Bottom 2-mL Assay Block (Costar, Corning Inc., Corning, New York, USA) was used with eight punched disks placed in each well. The disks were only added to the innermost wells (rows 3–10 out of 12) of the plates for more effective and cleaner transfer from one plate to another. Disk punches were taken from the leaves' darker print areas on the cards, which were presumed to contain the most concentrated residue.

Once disks had been collected from each plant, a series of washes were employed according to the manufacturer's protocol. First, 400 μL of FTA purification reagent was added to each well. The plate of samples was then covered with

Microseal 'F' foil seals (Bio-Rad, Hercules, California, USA) to reduce contamination, vortexed, and incubated at room temperature for 4 min before centrifugation at 6200 rpm for 2 min at 4°C. The supernatant was removed and the plates were then centrifuged a second time before the foil was removed. The FTA purification reagent wash was repeated once, followed by two low TE buffer washes. After completing the four washes, the remaining supernatant was removed. The samples were left to dry at room temperature for 20 min. The disks were then transferred to a new plate and centrifuged to separate them from as much remaining supernatant as possible. Any punches not in the new plate were manually transferred using cleaned forceps. Finally, 80 μL of TE was added to the plate containing the washed disks. The plate was centrifuged at 6200 rpm for 1 min and then incubated for 5 min at 95°C. The paper disks were left in the sample wells with the TE and eluted DNA.

CTAB-based extraction—For each replicate, a 1.0-cm²-sized piece of dried and shredded leaf tissue was added to a 2.0-mL tube containing 2.3-mm-diameter zirconia-silica beads and 1.0-mm-diameter glass beads (BioSpec, Bartlesville, Oklahoma, USA). The tissue was then macerated using the TissueLyser (QIAGEN) at 30 Hz for 30 s before 400 μL of CTAB, warmed to 65°C, was added. The Cactaceae samples produced a viscous gel, so 75 μL of 20 mg/mL Proteinase K and 500 μL of CTAB were added to further clean the sample. All samples were incubated overnight in a rotary incubator at 65°C and 150 rpm.

The resulting lysate solutions were centrifuged at 13,000 rpm for 5 min and 300 μ L of supernatant was transferred to each well of a 96-well AutoGen plate. DNA extraction was completed using the automated DNA isolation system AutoGenprep 965 (AutoGen) as outlined by the manufacturer's protocol, and the final pellets were resuspended in 80 μ L of TE buffer.

Quantification and quality assessment—To assess the quality of the extractions generated with the two techniques, 260/280 nm absorbance ratios and fluorometric determinations of DNA concentration were performed using a Synergy HT Microplate Reader (BioTek, Winooski, Vermont, USA). The Quant-iT Broad-Range dsDNA Assay Kit (Invitrogen, Waltham, Massachusetts, USA) was used with the Synergy Microplate Reader according to the manufacturer's protocol (Fig. 2). The DNA size range and degree of fragmentation were determined using gel electrophoresis (Fig. 3). The samples (10 μ L) were run alongside a HiLo DNA size standard (Minnesota Molecular, Minneapolis, Minnesota, USA) on a 1.5% SeaKem Agarose LE gel in 1× SB (Brody and Kern, 2004). The loading dye used to prepare these electrophoretic separations contained a 1:1000 dilution of Gel Red fluorescent nucleic acid stain (Biotium, Fremont, California, USA).

PCR amplification and DNA sequencing—To compare the amplification success of DNA extracted with the two methods, a portion of the low-copy nuclear gene At103, the nuclear ribosomal intergenic spacer (ITS), and the plastid ribulose-bisphosphate-dismutase large subunit (rbcL) were amplified. Each amplification reaction contained 2.5 μ L of FTA card–extracted sample or 2.5 μ L of 1:50 diluted CTAB-extracted sample. The At103 region was amplified with forward primer CTTCAAGCCMAAGTTCATCTTCTA and reverse primer

Table 1. Species sampled in this study. All specimens were collected in Camp Springs, Maryland, USA (GPS coordinates 38°50′40.4″N, 76°56′17.4″W). Triplicate vouchers were made for deposit at the National Museum of Natural History's Herbarium (US), the University of Illinois at Urbana-Champaign Plant Biology Herbarium (ILL), and the Chicago Botanic Garden (CHIC). Samples are organized alphabetically by family name.

Family	Genus/Species (Common name)	Voucher no.
Aquifoliaceae	Ilex glabra (L.) A. Gray (inkberry)	C. Siegel 11
Asclepiadaceae	Asclepias syriaca L. (common milkweed)	C. Siegel 10
Aspleniaceae	Asplenium platyneuron (L.) Britton, Sterns & Poggenb. (ebony spleenwort)	C. Siegel 6
Asteraceae	Ratibida pinnata (Vent.) Barnhart (pinnate prairie coneflower)	C. Siegel 8
Cactaceae	Opuntia cf. laevis J. M. Coult. (prickly pear)	C. Siegel 15
Cyperaceae	Carex lurida Wahlenb. (shallow sedge)	C. Siegel 3
Fabaceae	Albizia julibrissin Durazz. (silktree)	C. Siegel 7
Lamiaceae	Monarda fistulosa L. (wild bergamot)	C. Siegel 1
Magnoliaceae	Magnolia virginiana L. (sweetbay magnolia)	C. Siegel 12
Oxalidaceae	Oxalis dillenii Jacq. (wood sorrel)	C. Siegel 13
Pinaceae	Pinus virginiana Mill. (Virginia pine)	C. Siegel 5
Poaceae	Dichanthelium commutatum (Schult.) Gould (panicgrass)	C. Siegel 14
Simaroubaceae	Ailanthus altissima (Mill.) Swingle	C. Siegel 9
Typhaceae	Typha angustifolia L. (narrowleaf cattail)	C. Siegel 2
Vitaceae	Ampelopsis glandulosa (Wall.) Momiy. (Amur peppervine)	C. Siegel 4

http://www.bioone.org/loi/apps 2 of 7

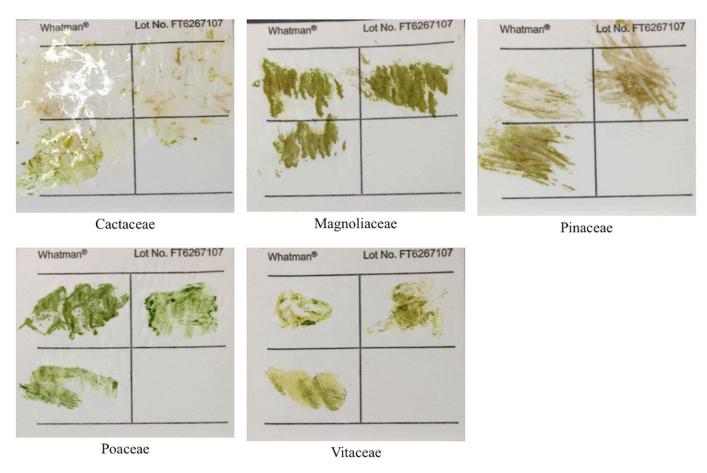


Fig. 1. Whatman FTA PlantSaver card samples, pre-extraction. The four quadrants of several FTA cards covering the phylogenetic range of the sampled species are shown. Three quadrants of each card were used to press leaf tissue. Eight hole-punches were obtained for each well of a 96-well plate.

TTGGCAATCCATTGAGGTACATNGTM (Li et al., 2008), and the ITS region was amplified with primers ITS5a (CCTTATCATTTAGAGGAAGGAG) and ITS4 (CAGGAGACTTGTACACGGTCCAG; Kress et al., 2005). For these primer pairs, the PCR program used was: 95°C for 4 min; four cycles at 95°C for 45 s, 57°C for 30 s, 72°C for 1.5 min; four cycles at 95°C for 45 s, 54°C for 30 s, 72°C for 1.5 min; 35 cycles at 95°C for 45 s, 52°C for 30 s; 72°C for 1.5 min; and a final extension at 72°C for 10 min. For *rbcL*, PCR was conducted with forward primer ATGTCACCACAAACAGAGACTAAAGC and reverse primer GTAAAATCAAGTCCACCRCG (Kress et al., 2005) using the following thermal cycler protocol: 95°C for 3 min; 34 cycles at 94°C for 30 s, 55°C for 45 s, 72°C for 2 min; and a final extension at 72°C for 5 min. Success of the PCR was assessed by gel electrophoresis.

Two successful replicates for each marker and species were chosen for sequencing. These products were enzyme-treated with ExoSAP-IT (Affymetrix, Santa Clara, California, USA) to remove unincorporated primer and deoxynucleotide triphosphates (dNTP) prior to being used as template for cycle sequencing with BigDye Terminator version 3.1 Ready Reaction Mix (Thermo Fisher Scientific, Grand Island, New York, USA). The resulting Sanger sequencing fragments were purified through Sephadex G-50 and analyzed on an ABI 3730xl automated capillary sequencer (Life Technologies, Carlsbad, California, USA). Chromatograms were analyzed using Geneious 9.0.5 (Biomatters, Auckland, New Zealand), and successful sequences were submitted to GenBank (Table 2).

RESULTS


Quantification and quality assessment—Spectrophotometric estimations of nucleic acid concentrations of CTAB-extracted samples were considerably higher than those of FTA


card–extracted samples (Fig. 2A). However, there was not a substantial difference in average detected concentration among the 15 sampled species. Although the CTAB-extracted samples were much more concentrated, the CTAB sample concentration values were more variable than those of the FTA card–extracted samples. FTA card samples had an average standard error for concentration of 7.81 ng/ μ L while that of the CTAB-extracted samples was 37.96 ng/ μ L. The average 260/280 absorbance ratio was closer to 1.8 for CTAB-extracted samples compared to those for FTA card–extracted samples (Fig. 2B). The FTA card–extracted, amplified samples also had 260/280 ratios that were indicative of lower quality.

Quant-iT fluorometry showed that the CTAB-extracted samples contained large, but highly variable concentrations of DNA, while the DNA extracted with the FTA cards was relatively undetectable with this method (Fig. 2C). Like the spectrophotometric measurements, there were no substantial differences in average concentration among the 15 sampled species. Because spectrophotometry does not discriminate between RNA and DNA as effectively as the fluorometric method, the greater quantities of DNA estimated by the spectral absorbance values (Fig. 2A vs. Fig. 2C) might be an artifact of the copurified RNA with the DNA in the samples.

Gel separations of the samples indicate that all of these total genomic DNA extracts contained unfragmented, high-molecular-weight DNA, except for the Typhaceae samples extracted with the FTA card method (Fig. 3).

http://www.bioone.org/loi/apps 3 of 7

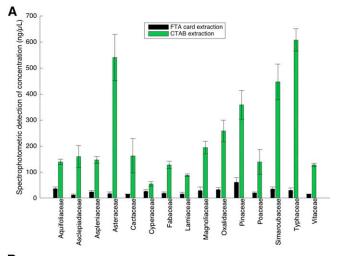

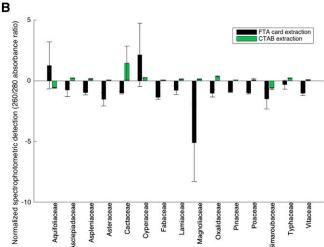

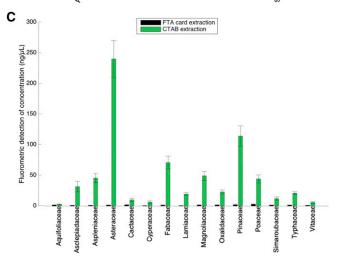


Fig. 2. DNA concentration and quality in plant extractions. (A) Concentration ($ng/\mu L$) \pm 1 SE as detected by the Epoch Spectrophotometer System determined for seven replicates. Nucleic acid concentration was indicated by absorbance measured at 260 nm. Measurements did not discriminate between RNA and DNA concentrations. (B) 260/280 absorbance ratio as detected by the Epoch Spectrophotometer System. Absorbances measured at 260 nm and 280 nm with suboptimal 1.8 absorbance ratios indicated the presence of contaminants. Detected values have been normalized to 0 by subtracting 1.8 from each measured ratio value. (C) Concentration

PCR amplification and DNA sequencing—CTAB and FTA card samples for all families were generally successful in amplification of At103, ITS, and rbcL (Table 2). Despite the differences seen in the appearance of each plant print on the FTA cards (Fig. 1), all but the conifer *Pinus virginiana* could be successfully amplified with at least one primer pair (Table 2). The rbcL amplicons were the most amenable to cycle sequencing, resulting in highquality chromatograms for at least one CTAB and/or FTA card sample of each species (Table 2). The At103 marker most readily amplified and sequenced for families Aquifoliaceae, Aspleniaceae, Cactaceae, Fabaceae, Lamiaceae, Oxalidaceae, Poaceae, Simaroubaceae, Typhaceae, and Vitaceae. The most successful families for amplification and sequencing with the ITS5a/ITS4 primers were Aquifoliaceae, Asclepiadaceae, Asteraceae, Cactaceae, Cyperaceae, Poaceae, Simaroubaceae, and Vitaceae. These results indicated that, even if the leaf did not create a dark chlorophyll print, amplifiable DNA was captured by the FTA card.

DISCUSSION

The FTA card method could have great utility in the study of non-agricultural plant phylogenetics and population genetics, as it addresses some shortcomings of the CTAB-based technique, including facility of collection and transportation. During field collection, kilogram quantities of silica would not need to be transported to preserve specimens for CTAB-based extraction methods, allowing more specimens to be collected. While collection permits are usually required for sampling plants, a U.S. Department of Agriculture Animal and Plant Health Inspection Service (USDA-APHIS) permit is not necessary for acquiring the actual plant tissue embedded on an FTA card, because the material cannot be propagated, as would whole-rooted plants, rhizomes, and seeds (V. Funk, personal communication). Additionally, because the FTA card embedding procedure may destroy RNA from viruses (Kraus et al., 2011), biosafety issues may not arise when transporting plant tissues between states and between countries. The FTA cards are more compactly stored and do not require refrigeration. Finally, the FTA card-extraction method requires less laboratory expertise and fewer hazardous chemicals such as CTAB, chloroform, and phenol (Suzuki et al., 2006; Marques et al., 2010; Chandrashekara et al., 2012).

The FTA card and CTAB extraction methods both exhibited varying levels of success. CTAB-extracted samples contained higher concentrations of DNA as estimated by their A260/A280 absorbance ratios. Between replicates of the same species there was greater variation in concentration among those extracted with the CTAB-based vs. the FTA card methods. This inconsistent quantity of DNA recovered from the CTAB procedure may have been an artifact of the AutoGen instrument, as has been previously reported (Mulcahy et al., 2016).

The increased concentrations of DNA detected in the CTAB-extracted samples is a major consideration in determining the overall utility of the two methods. The differences in DNA concentration for CTAB- and FTA card-extracted samples could have been due to the differences in amount of leaf tissue originally used for extraction. Approximately 1 cm² of leaf tissue was used for each CTAB-extracted sample, while for FTA

 $(ng/\mu L)$ as detected by the Quant-iT Fluorometer. Fluorescence-based dyes were bound to the DNA in each sample. This technique does discriminate between RNA and DNA.

http://www.bioone.org/loi/apps 4 of 7

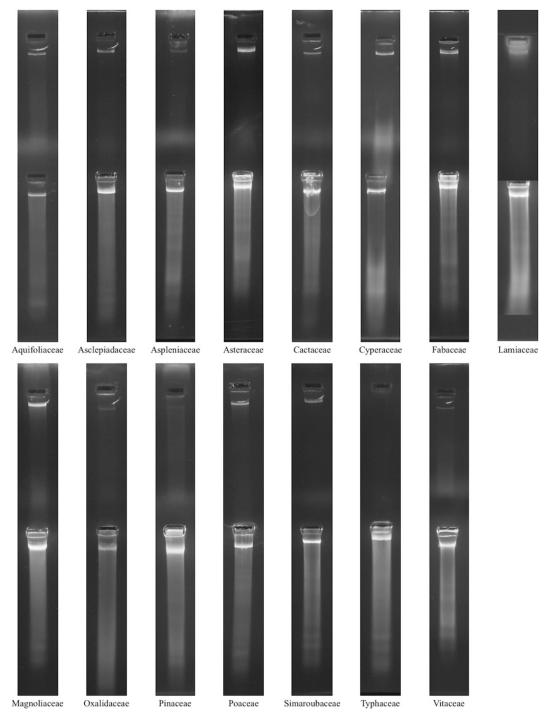


Fig. 3. Total genomic DNA extract gel electrophoresis results. The FTA card-extracted DNA is seen in the top row while the CTAB-extracted DNA is seen in the bottom row. All seven replicates of each plant extraction method were run on the gel, but only one replicate per species is shown.

card-extracted samples the amount of DNA produced depended on the concentration and amount of fluid transferred from the leaves to the card.

The total genomic DNA gels (Fig. 3) indicated that the FTA card–extracted samples contained a less-fragmented DNA when compared to the CTAB-extracted samples. In contrast to the FTA card–extracted DNA, it was observed that RNA and/or degraded DNA could easily be detected in gel separations of

the CTAB-extracted DNA. The amount of PCR product generated was approximately correlated with the concentration of the DNA extract used. Although many of the amplification success rates were similar between the two extraction methods, the CTAB-extracted samples on average amplified more consistently (Table 2). Only the Cactaceae, Cyperaceae, and Simaroubaceae species were amplified more successfully with the FTA card-extracted samples. In addition to the conifer and

http://www.bioone.org/loi/apps 5 of 7

Table 2. Results of DNA amplifications and sequencing.^a

			Extracted	Extracted using CTAB method	ethod		Extracte	Extracted using FTA method	thod
Family	Marker	No band	>1 band	1 band	GenBank accession no. ^b	No band	>1 band	1 band	GenBank accession no. ^b
Aquifoliaceae	At103	0	7	0	No	0	7	0	KX394243
	SLI	0	0	_	KX352749	0	0	7	KX352750
	rbcL	4 -	0 (m t	1;	0 (7	0 0	KX702293
Asciepiadaceae	AHUS	4 <	7 1	٦ ٥	NO 1NO 1NO	~ 0	7 0	0 0	NO 100 100 100 100 100 100 100 100 100 10
	rbcI	0 0	~ C	7 0	# 1255XX	0 0	· C	7 0	CX702387
Asnleniaceae	At 103	9	0 0		KX 394240	0 1	0 0	· C	(877)07781 —
rispiranteed	SLI	0 0	v		OF TON	-	-	v	Z
	rbcL	0	0	1 1-	KX702288	2	0	າທ	KX702289
Asteraceae	At103	2	ν.	0	oZ	8	-	n	No.
	SLI	0	0	7	KX352753	0	0	7	KX352754
	rbcL	0	0	7	KX702303	0	5	7	KX702304
Cactaceae	At103	0	7	0	KX394245	0	0	7	KX394246
	ITS	0	7	0	KX352751	1	9	0	KX352752
	rbcL	1	0	9	KX702298	0	7	0	KX702299
Cyperaceae	At103	7	0	0		9	0	-	No
4	ITS	1	1	w	No	0	0	7	KX352746
	rbcL	0	9	1	KX702290	0	2	w	KX702291
Fabaceae	At103	0	7	0	KX394236	0	7	0	KX394237
	ITS	0	9	1	No	0	7	0	No
	rbcL	0	2	w	KX702282	0	_	9	KX702283
Lamiaceae	At103	0	0	7	KX394244	9	1	0	No
	SLI	1	4	7	No	7	0	0	
	rbcL	0	7	0	KX702296	4	0	3	KX702297
Magnoliaceae	At103	7	0	0	ı	S	2	0	No
	SLI	0	7	0		0	7	0	
	rbcL	0	5	7	KX702295	0	7	0	KX702294
Oxalidaceae	At103	0	7	0	KX394247	Т	5	_	KX394248
	SLI	0	0	7	No	0		9	No
	rbcL	0	4	က	KX702300		9	0	KX702301
Pinaceae	At103	0	7	0	No	7	0	0	ı
	SLI	0	5	7	No	7	0	0	1
	rbcL	0	0	7	KX702302	7	0	0	
Poaceae	At103	0	7	0	KX394241	3	3	_	KX394242
	SLI	0	7	0	KX352747	0	7	0	KX352748
	rbcL	0	4	8	No	2	5	0	KX702292
Simaroubaceae	At103	3	4	0	KX394234	0	0	7	KX394235
	SLI	2	3	2	No	0	7	0	KX352741
	rbcL	3	0	4	KX702280	0	7	0	KX702281
Typhaceae	At103	0	0	7	KX394249	3	0	4	KX394250
	SLI	0	9	1	So	0	7	0	No
	rbcL	0	4	က	KX702305	2	_	4	No
Vitaceae	At103	0	0	7	KX394238	4	0	က	KX394239
	SLI	_	1	w	KX352742	0	2	7	KX352743
	rbcL	0	0	7	KX702284	4	ec	0	KX702285
a Domilion one line	- Inchohoticolis	a Formiliae are listed alphabationally with the mumber of renaise		of of a total of	otes (mit of a total of seven) that shows and and	multiple bonde or	band one	The minutes of wallacted that was	The first work of a color

*Families are listed alphabetically with the number of replicates (out of a total of seven) that showed no bands, multiple bands, or one band. The number of replicates that were successfully amplified with a single band are provided in boldface. Plant samples were amplified with three different primer pairs: ITS, *rbcL*, and *At103*.

*b — = PCR was not successful and, as a result, sequencing was not performed; No = sequencing was performed on an amplicon, but was not successful. GenBank accession numbers are provided for amplicons for which PCR and sequencing were successful.

fern species, the families of the species that successfully amplified and sequenced were distributed evenly throughout the Angiosperm Phylogeny Group (APG) III phylogeny for flowering plants (Angiosperm Phylogeny Group, 2009; Tables 1, 2).

In summary, this study provided a more comprehensive understanding of the strengths and weaknesses of using Whatman FTA cards for non-agricultural species. Although the FTA cardextracted samples had a substantially lower concentration of DNA, the method is a good alternative for field collection if well-optimized primer pairs are being used that can amplify low concentrations of template DNA.

The previous success of FTA cards with large sample sizes in agricultural studies indicates that it could be a suitable method for studies in plant population genetics and conservation biology that require the analysis of hundreds to thousands of samples (Adugna et al., 2011; Chandrashekara et al., 2012). The FTA card method may work better for some families than for others. For example, we believe the FTA card method may work well for those families and genera phylogenetically related to and/or physiologically similar to the Cactaceae and Pinaceae.

Finally, while FTA cards have been used previously for extracting insect nucleic acids for RNA-Seq next-generation sequencing studies (Grubaugh et al., 2015), future work should focus on whether plant DNA and/or RNA extracted with the FTA cards could be used for some next-generation sequencing techniques (i.e., RNA-Seq, and possibly RAD-Seq, Genotyping by Sequencing [GBS], and targeted enrichment; Zimmer and Wen, 2015). By increasing the number and size of punches produced from the cards, FTA card–extracted samples might provide enough material suitable for the latter DNA-based next-generation sequencing methods.

LITERATURE CITED

- ADUGNA, A., P. M. SWEENEY, AND A. A. SNOW. 2011. Optimization of a high throughput, cost effective, and all-stage DNA extraction protocol for sorghum (*Sorghum bicolor*). *Journal of Agricultural Science and Technology* 5: 243–250.
- Angiosperm Phylogeny Group. 2009. An update of the angiosperm phylogeny group classification for the orders and families of flowering plants: APG III. *Botanical Journal of the Linnean Society* 161: 105–121.

- Brody, J. R., and S. E. Kern. 2004. Sodium boric acid: A Tris-free, cooler conductive medium for DNA electrophoresis. *BioTechniques* 36: 214–216.
- Bujang, N. S., N. A. Harrison, and N. Su. 2011. An improved method for extraction and purification of termite Endo- β -1,4-Glucanase from FTA® Cards. *Florida Entomologist* 94: 356–358.
- Bustin, S. A., and T. Nolan. 2004. Pitfalls of quantitative real-time reverse-transcription polymerase chain reaction. *Journal of Biomolecular Techniques* 15: 155–166.
- CHANDRASHEKARA, K. N., M. K. PRASANNAKUMAR, M. DEEPA, AND A. VANI. 2012. A rapid, sensitive and reliable method for detecting *Ralstonia solanacearum* using FTA (Whatman) Card. *Journal of Plant Pathology* 94: 219–221.
- DOYLE, J. J., AND J. L. DOYLE. 1987. A rapid DNA isolation procedure for small quantities of fresh leaf tissue. *Phytochemical Bulletin* 19: 11–15.
- GRUBAUGH, N. D., S. SHARMA, B. J. KRAJACICH, L. S. FAKOLI, F. BOLAY, J. W. DICLARO, W. E. JOHNSON, ET AL. 2015. Xenosurveillance: A novel mosquito-based approach for examining the human-pathogen land-scape. PLoS Neglected Tropical Diseases 9: e0003628.
- Kraus, R. H. S., P. van Hooft, J. Waldenström, N. Latorre-Magalef, R. C. Ydenberg, and H. H. T. Prins. 2011. Avian influenza surveillance with FTA cards: Field methods, biosafety, and transportation issues solved. *Journal of Visualized Experiments* 54: e2832.
- KRESS, W. J., K. J. WURDACK, E. A. ZIMMER, L. A. WEIGT, AND D. A. JANZEN. 2005. Use of DNA barcodes to identify flowering plants. *Proceedings* of the National Academy of Sciences, USA 102: 8369–8374.
- LI, M., J. WUNDER, G. BISSOLI, E. SCARPONI, S. GAZZANI, E. BARBARO, H. SAEDLER, AND C. VAROTTO. 2008. Development of COS genes as universally amplifiable markers for phylogenetic reconstructions of closely related plant species. *Cladistics* 24: 727–745.
- MARQUES, A. P., L. ZÉ-ZÉ, M. V. SAN-ROMÃO, AND R. TENREIRO. 2010. A novel molecular method for identification of *Oenococcus oeni* and its specific detection in wine. *International Journal of Food Microbiology* 142: 251–255.
- Mulcahy, D. G., K. S. MacDonald, S. G. Brady, C. Meyer, K. B. Barker, and J. Coddington. 2016. Greater than *X* kb: A quantitative assessment of preservation conditions on genomic DNA quality, and a proposed standard for genome-quality DNA. *PeerJ Preprints* 4: e2202v1. https://doi.org/10.7287/peerj.preprints.2202v1.
- SUZUKI, S., H. TAKETANI, K. KUSUMOTO, AND Y. KASHIWAGI. 2006. High-throughput genotyping of filamentous fungus Aspergillus oryzae based on colony direct polymerase chain reaction. Journal of Bioscience and Bioengineering 102: 572–574.
- ZIMMER, E. A., AND J. WEN. 2015. Using nuclear gene data for plant phylogenetics: Progress and prospects II: Next-gen approaches. *Journal of Systematics and Evolution* 53: 371–379.

http://www.bioone.org/loi/apps 7 of 7